



Safety IV Cannula

T-1046

- Passive fully automatic protection.
- Needle sharp tip is enclosed inside a plastic cage after pulling out.
- The plastic cage is designed as such that once needle sharp tip is fully enclosed inside it can't come out of the cage hence prevent.
- Healthcare professionals from the risk of accident due to needle stick injury of injection.
- The size of Safety I.V. Cannula is very similar to conventional I.V. Cannula and the technique of cannulation is similar to cannulation technique of conventional I.V. Cannula.
Sterile / Disposable / Individually Packed.

Size	Color	Chtheter	Size	Length	Flow Rate
14 G	Orange	1.7	2.1	45	270
16 G	Grey	1.3	1.7	45	180
17 G	White	1.1	1.5	45	125
18 G	Green	0.9	1.3	45	80
20 G	Pink	0.8	1.1	32	54
22 G	Blue	0.6	0.9	25	33
24 G	Yellow	0.5	0.7	19	20
26 G	Violet	0.4	0.6	19	15

IV Cannula Fixator

T-1047





- Used for Fixation of I.V Cannula / Catheter.
- Made from Thin Polyurethane film and available in two variants -Transparent and Non-woven.
- Breathable and water proof, with reinforced notches that reduce mechanical stress from heavy catheters.
- Thin film ensures easy fixation of devices with contoured surfaces.
Sterile / Disposable / Individually Packed.

Scalp Vein Set

T-1050 [With or Without Luer Lock]

- Butterfly set for long term infusion.
- Scalp vein set is designed to provide rapid venous access with greater patient comfort during infusion.
- Short beveled siliconized needle facilitates atraumatic cannulation.
- Thin wall needle provides better flow rate per gauge, as in the same cannula more circumference is available for better flow of fluid.
- Butterfly shaped wings facilitate easy handling and attachment with the skin.
- The proximal end of the set is provided with flexible female luer fitting.
- Butterflies are color coded for instant identification of needle size.
- Butterfly is connected to soft non toxic, non-irritant, medical grade tube, which does not kink or coil set.
- Available with DEHP Free Material.
Sterile / Disposable / Individually Packed.

Color Code	Pink	Brown	Yellow	Green	Black	Sky Blue	Purple	Orange
Size G	18	19	20	21	22	23	24	25

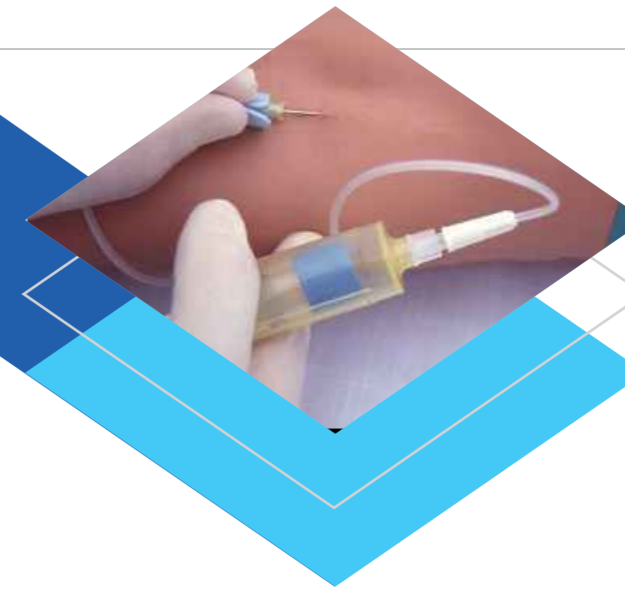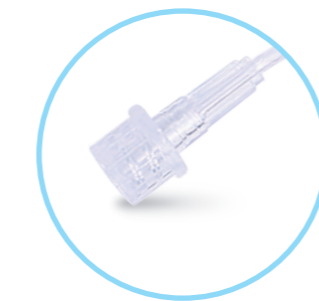


Measured Volume Burette Set

T-1051 to T-1054

- Super smooth kink resistance tubing ensures uniform flow rate.
- Flexible, Transparent & soft cylindrical measured volume chamber with sharp and easy piercing Spike.
- Smooth Roller Clamp Facilitates easy, safe control and adjustment of fluid rates.
- Specially designed floating auto shut off valve acts as floating indicator and automatically shut off the drain path when the chamber gets empty to prevent air trapping in the fluid line.
- Separate plug for extra medication and continuous change over.
Sterile / Disposable / Individually Packed.

OPTIONS AVAILABLE

- With or Without Airvent.
- With or Without Needle.
- With or Without "Y" Connection.
- With or Without Luer Lock.
- Latex or Latex Free.
- Available with DEHP Free Material.
- Tube Length : 150cms, 180cms, 200cms.
- Size : 150ml., 110ml.

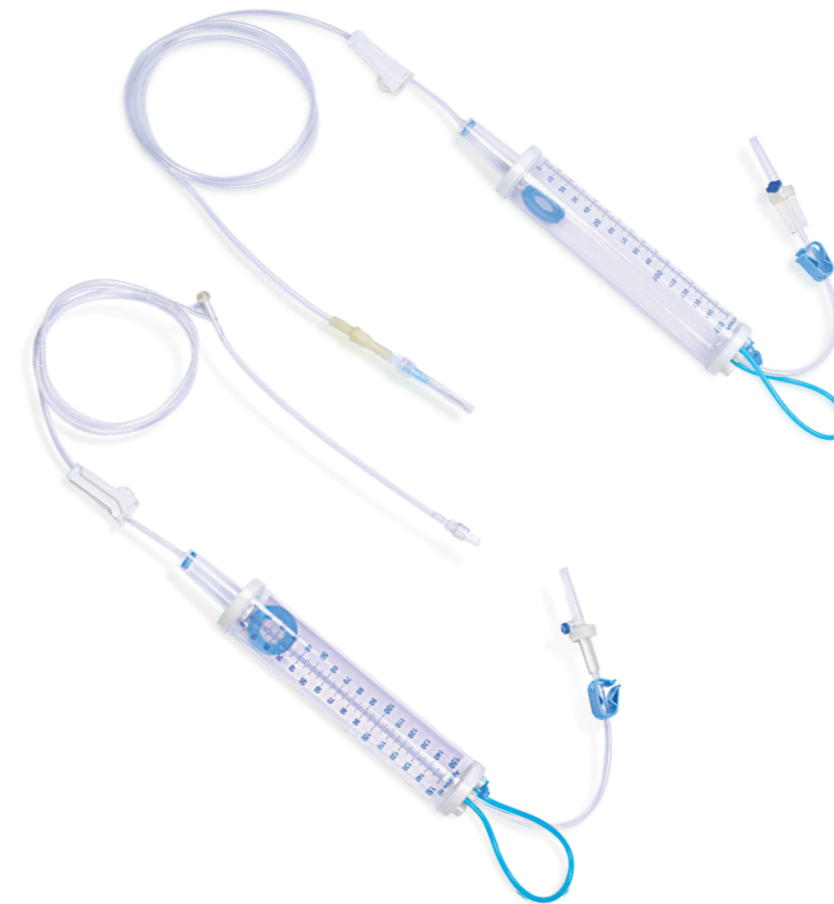

Extension Tube

T-1056

- Super smooth kink resistance tubing ensures uniform flow rate.
- Male luer lock connector at one end and female luer lock connector at other end.
- Available with DEHP Free Material.
Sterile / Disposable / Individually Packed

Tube Length CM
10 25 50 80 100 150 200 500

OPTIONS AVAILABLE

- With or Without "Y" connection.
- With or Without Clamp.
- Customized Inner & Outer Diameter.

