

Endotracheal Tube

Cuffed & Plain
A - 3012




- Used to maintain a patent airway and mechanically ventilate patients.
 - High volume, low pressure cuff, maintains seal throughout ventilator cycle.
 - X-ray detective line available.
 - Thermosensitive PVC allows the tube to soften at body temperature to conform to airway anatomy.
 - Clear PVC tube with visible markings for easy observation.
 - Murphy Eye smoothly formed to allows ventilation in the event of obstruction of the end of the tube during intubation.
 - A full range of standard sized to meet most of your requirements.
 - Universal 15 mm Connector at proximal end.
 - Available with DEHP free Material.
- Sterile / Disposable / Individually Packed.*

	Size MM
CUFFED :	4.0, 4.5, 5.0, 5.5, 6.0, 6.5, 7.0, 7.5, 8.0, 8.5, 9.0, 9.5, 10.00
PLAIN :	2.0, 2.5, 3.0, 3.5, 4.0, 4.5, 5.0, 5.5, 6.0, 6.5, 7.0, 7.5, 8.0, 8.5, 9.0, 9.5, 10.00

Reinforced

Cuffed & Plain
A - 3013

- Used to maintain a patent airway and mechanically ventilate patients.
 - Uniform spiral metal-reinforced tube wall to prevent the tube from kinking.
 - Reinforcement supports tube position and verification during X-ray.
 - Accurate and safe positioning ensured by graduation marks and tip-to-tip radiopaque line.
 - High volume / low pressure cuff ensures efficient low pressure cuff seal (For cuffed tubes).
 - Universal 15 mm connector at proximal end.
 - Available with DEHP free Material.
- Sterile / Disposable / Individually Packed.*

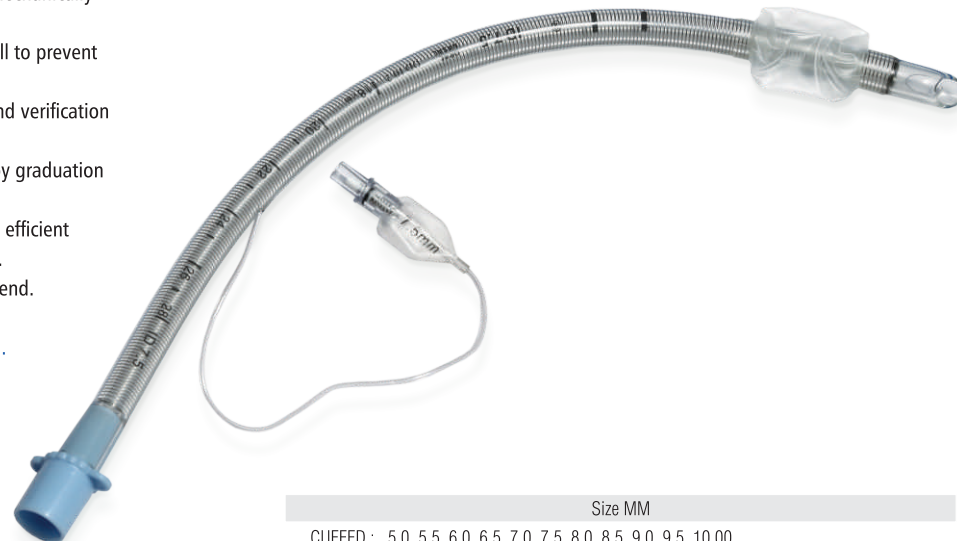

	Size MM
CUFFED :	5.0, 5.5, 6.0, 6.5, 7.0, 7.5, 8.0, 8.5, 9.0, 9.5, 10.00
PLAIN :	3.0, 3.5, 4.0, 4.5, 5.0, 5.5, 6.0, 6.5, 7.0, 7.5, 8.0, 8.5, 9.0, 9.5, 10.00



L. P. Spinal Needle

A - 3014

Quincke Point or Pencil Point
With or Without Introducer

- Spinal needle provided exceptional control when penetrating the dura.
 - Suitable for spinal anaesthesia.
 - Aggressive anesthetic distribution upon injection.
 - High flow rate enables faster cerebro spinal (CSF) flashback.
 - Clear polycarbonate hub offer easy visualization of CSF or blood.
 - Available with Customized Needle Length.
- Sterile / Disposable / Individually Packed.*



Color Code	Pink	White	Yellow	Green	Black	Light Blue	Violet	Orange	Brown	Grey
Size G	18	19	20	21	22	23	24	25	26	27

Nebulizer Mask

A - 3015 - 3016

- Designed to be used with any standard Nebulizer kit.
 - Super smooth kink resistance tubing ensures uniform flowrate.
 - Smooth and feathered edge of face mask for patient comfort and reducing irritation points.
 - Adjustable elastic strip and integrated nose clip for proper positioning of mask.
 - The Nebulizer Kit just plugs into the standard fitting located at the bottom of the Mask.
 - Masks are designed to deliver more medication directly to the patient due to the direct flow design.
 - Proximal end of tube is connected with funnel shape connector for easy connection with oxygen source.
 - Star lumen main tube to avoid accidental blockage.
 - Tube Length : 210cm.
 - Size : Adult, Paediatric & Neonatal.
 - Available with DEHP free Material.
- Sterile / Disposable / Individually Packed.*

