

### T-PC. For Nebulizer

A - 3024



- Super smooth kink resistance tubing ensures uniform flowrate.
- Smooth and feathered edge of face mask for patient comfort and reducing irritation points.
- Quick to assemble - Facilitates fast nebulisation
- Can be used with Mask, iGel, LMA and ET tube
- Rotatable sections - Can be used on patients in sitting or supine position.
- Drug dosages and oxygen flow rates remain as per guidelines.
- Proximal end of tube is connected with funnel shape connector for easy connection with oxygen source.
- Star lumen main tube to avoid accidental blockage.
- Tube Length : 210cm.
- Size : Adult, Paediatric & Neonatal.
- Available with DEHP free Material.  
Sterile / Disposable / Individually Packed.



### Catheter Mount

A - 3025 [ Standard / Expandable ]


- Catheter Mount connects patient ET to breathing circuit. This minimizes the jerk of breathing to ET or TT reducing trauma to trachea.
  - The double swivel makes no dragging on breathing circuit & ET or TT.
  - Light weight, Double swivel rotates 360° on both axis, Minimum compliance, Low dead.
  - Compatible with all types of breathing and ventilator circuits.
  - Provided with 22 mm standard female connector on both the ends with inner diameter of 15 mm.
  - Collapsible corrugated tubing being inert to all anesthetic gases and reagents.
  - An elbow port is provided for easy suctioning without causing disturbance to breathing circuits.
- Sterile / Disposable / Individually Packed.



### Bain Circuits

A - 3026

- Bain Circuit offered comprises co-axial modification of basic T-piece system which has been developed for facilitating scavenging of waste anesthetic gases. As a tube carrying fresh gas, it travels inside outer reservoir tube to endotracheal tube connector. The process includes patient inspiring fresh gas from the outer reservoir tube and expiring into reservoir tube.
  - Adult : It is based on Mapleson D system, provided with 22 mm corrugated tube, expiratory valve and 2 litre antistatic bag.
  - Pediatric : It is based on Mapleson F system, provided with 15 & 10 mm corrugated tube with a 0.5 ltr antistatic Rebreathing Bag and a bag bleed valve.
- Sterile / Disposable / Individually Packed.

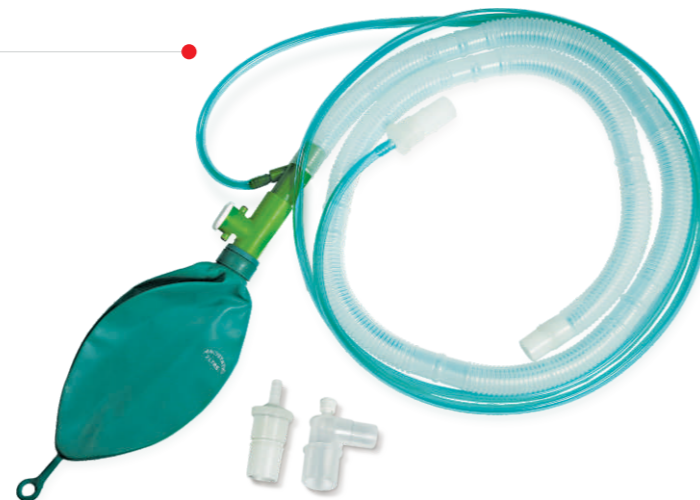


### Ventilator Circuits

A - 3027


[ Plain, Expandable, Single Water Tap, Double Water Tap ]

- Ventilator Circuit is the gas pathway connecting patients lung to the mechanical ventilator. It is essential that Circuit should be highest quality in order to ventilate patient correctly, efficiently and safely. We offer wide range of circuit with ISO standard connectors to suit most of the ventilator and humidifier.
- The corrugated tubes have a good flexibility and bending resistance, not meander, not broken, that mechanical ventilation in the process smooth gas supply.
- Breathing Circuit uses standard connector that is simple convenience in clinic.
- Size : Adult, Paediatric & Neonatal.  
Sterile / Disposable / Individually Packed.

Single Water Tap

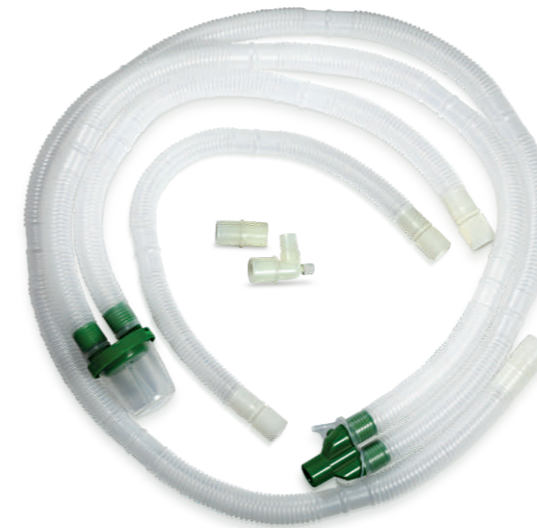

Plain

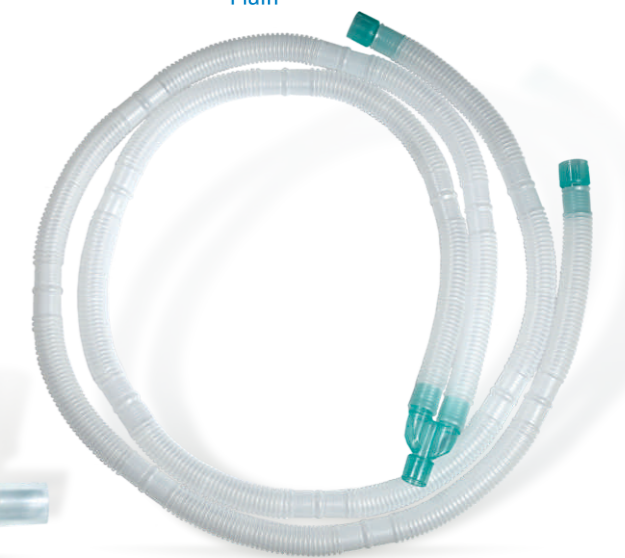

Options Available

- With or Without Port.
- With or Without "Y" - Piece.

Expandable

