


### Levin's Tube

G - 5003

- Suitable for gastro intestinal feeding & aspiration.
- Super smooth kink resistance tubing ensures uniform flowrate.
- Atraumatic, soft rounded, closed tip with four lateral eyes for efficient aspiration & administration.
- Tube is marked at 45,55,65 and 75 cm from the tip for accurate placement.
- Proximal end is fitted with universal funnel shaped connector for extension.
- Color coded connectors for easy identification of size as per standards.
- Available with DEHP Free Material.
- Available with universal luer Connector.
- Tube Length: 125 cm.

Sterile / Disposable / Individual Packed.

Color Code	Light Green	Blue	Black	White	Green	Orange	Red	Yellow	Violet	Light Blue
Size FG	6	8	10	12	14	16	18	20	22	24

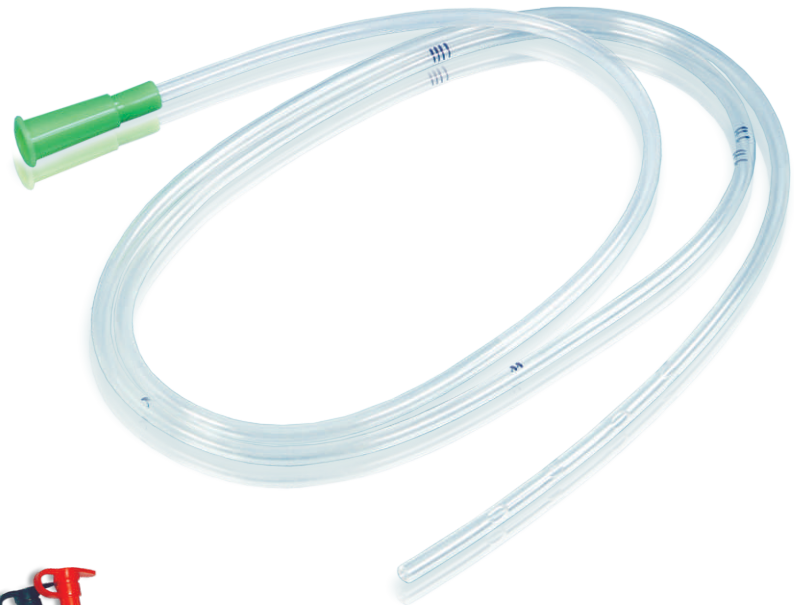


### Infant Feeding Tube

G - 5004

- For neonatal and paediatric nutritional feeding.
- Super smooth kink resistance tubing ensures uniform flowrate.
- Atraumatic, soft rounded, closed tip with two lateral eyes for efficient aspiration & administration.
- Proximal end is fitted with female luer mount for easy connection to feeding funnel or syringe.
- Color coded connectors for easy identification of size as per standards.
- Tube is marked at 20cm away from the distal tip.
- X-Ray opaque line provided through out the length of Tube.
- Available with graduation mark after every cm.
- Available with DEHP Free Material.
- Tube Length : 52cm.

Sterile / Disposable / Individual Packed.

Color Code	Grey	Green	White	Blue	Orange	Black
Size FG	5	6	7	8	9	10

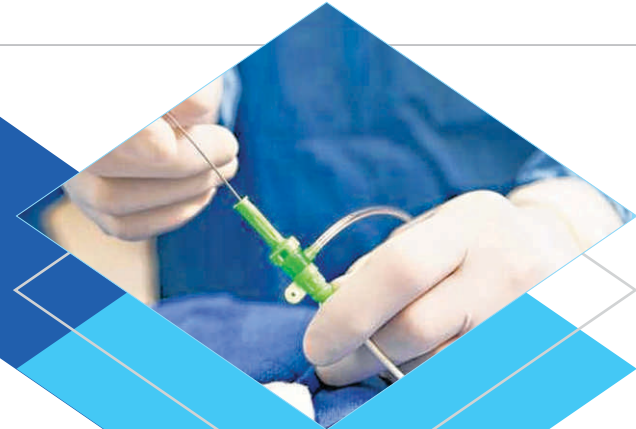

### Colostomy Kit

G - 5006

- The colostomy kit gives perfect relief and comfort during prolonged use.
- It is light weight, flexible and the belting arrangement fits the body well.
- 15mm 'U' form padding is provided between the sleeve and main body for better comfort.

Non-Sterile / Individually packed.



### Yankauer Suction Set

G - 5007 [ Vaccu Suck Set ]

- Yankaur Suction Set is suitable for convenient removal of secretion, blood and debris etc per operatively.
- Complete set is provided with Yankaur Suction Tip mounted on two meter long ribbed tube provided with Universal connector at both ends.
- Universal connectors are moulded from soft PVC, so as to accommodate the Suction tip on one end and all type of connectors of suction apparatus at the other end.
- Ribbed tube is kink resistant and has the strength to withstand the suction.

Sterile / Disposable / Individually Packed.



### Umbilical Catheter

G - 5008

- Designed for intermittent or continuous access to the umbilical artery or vein of newly born or premature infants.
- Super smooth kink resistance tubing ensures uniform flowrate.
- Smooth round tip helps easy cannulation and silk finish surface further facilitates the passage of catheter into the vein.
- Open distal end without lateral eyes, eliminates the chance of clot formation in the blind spaces.
- X-Ray opaque line provided through out the length of Tube.
- Clear scale marking at every cm helps to ascertain depth of insertion.
- Proximal end is fitted with female luer mount for convenient connection to I.V. therapy equipment.

Sterile / Disposable / Individually packed.

Color Code	Grey	Green	Blue
Size FG	5	6	8

